



Turn on

To power the CES device, the user presses and holds the power button until the power LED switches on, the battery is briefly displayed on the numbers labelled 1 to 8.

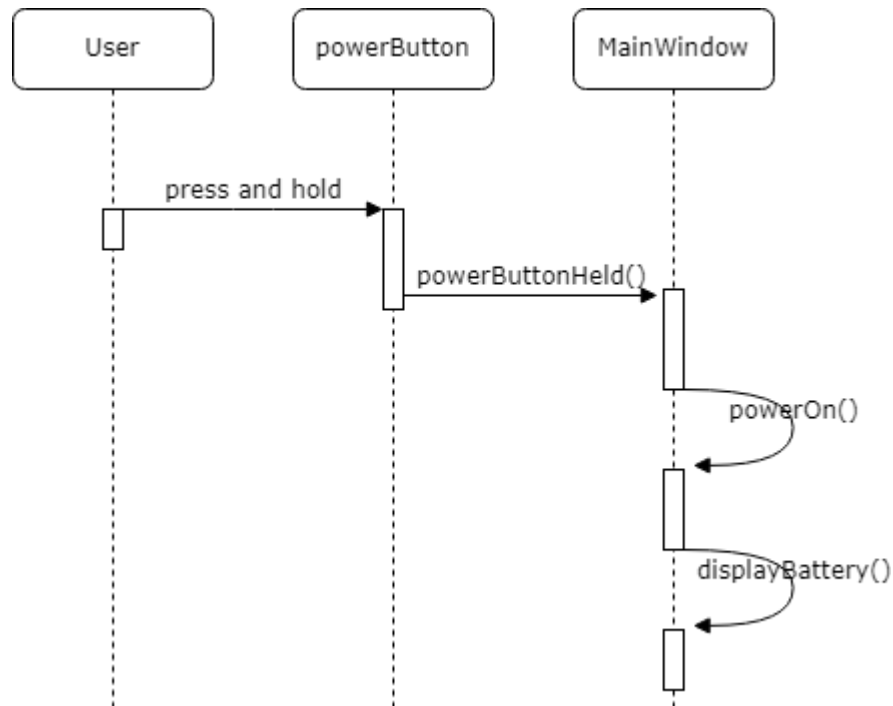


Turn off

If there are no active sessions and the power button is held while the device is switched on, the device turns off.

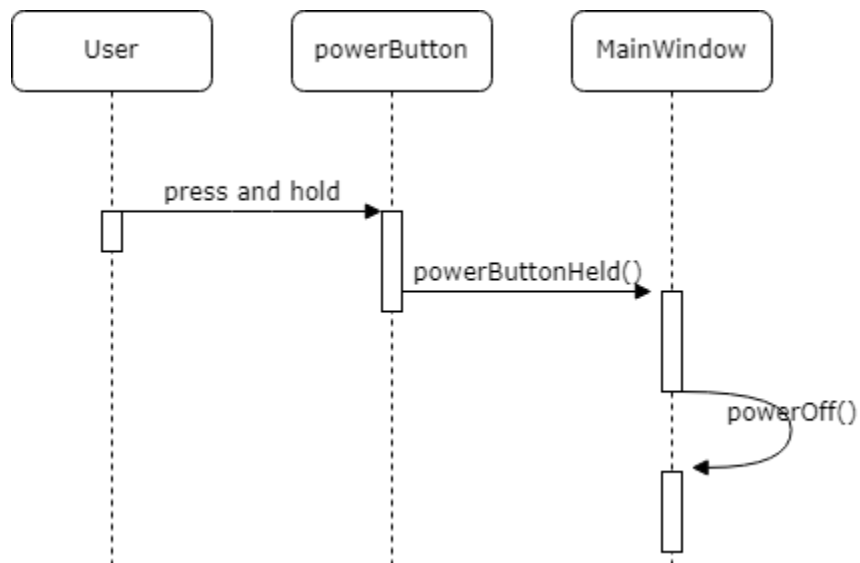


Selecting Sessions

To select session, the user presses the power button to cycle between the groups, then uses the INT buttons to select between the different session types

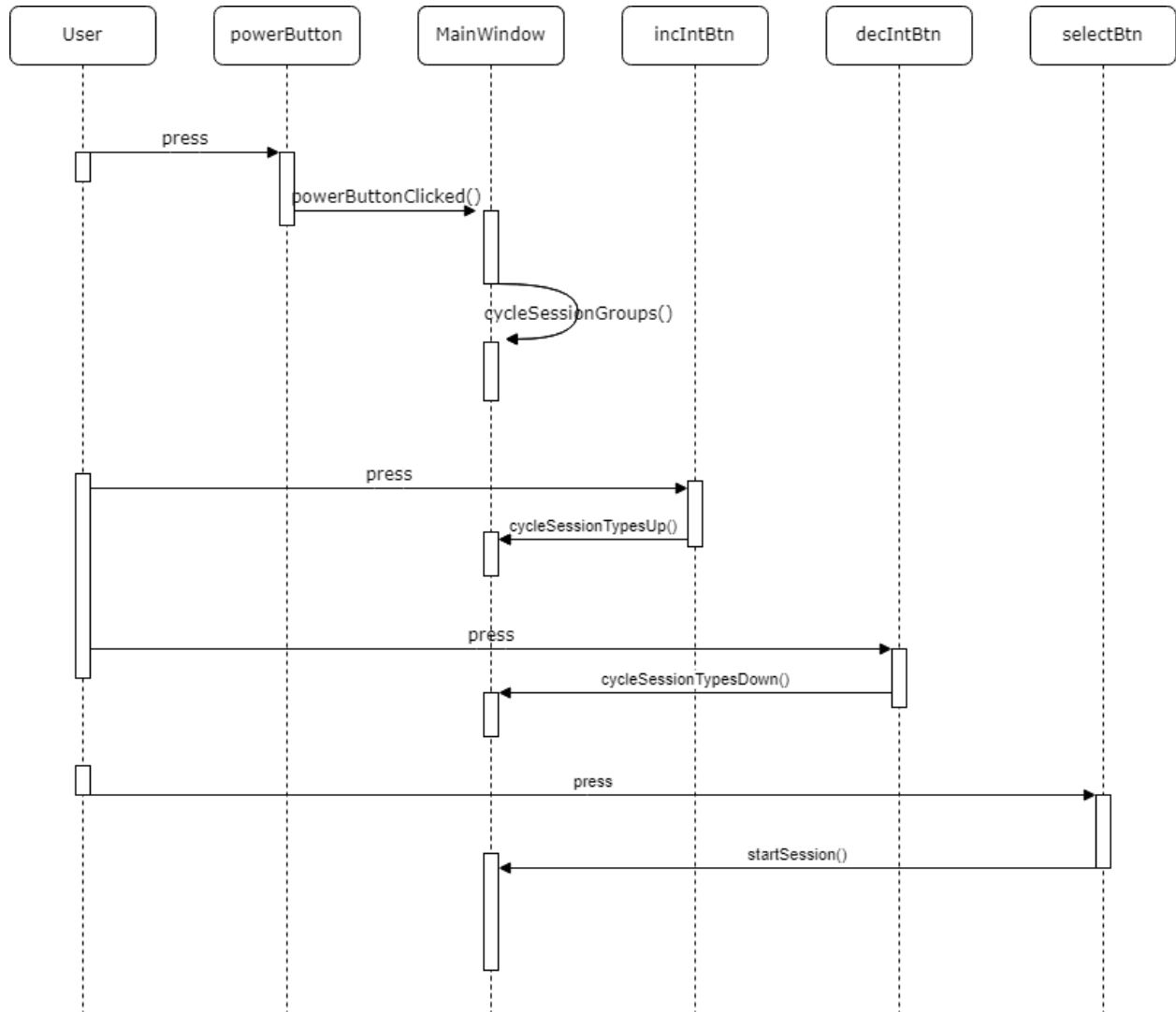


Connection Test

To begin a connection test, a session must be started by pressing the start session button. The device will then flash the CES Mode indicator and tDCS number corresponding to the selected session. After the test is done, the session will start.

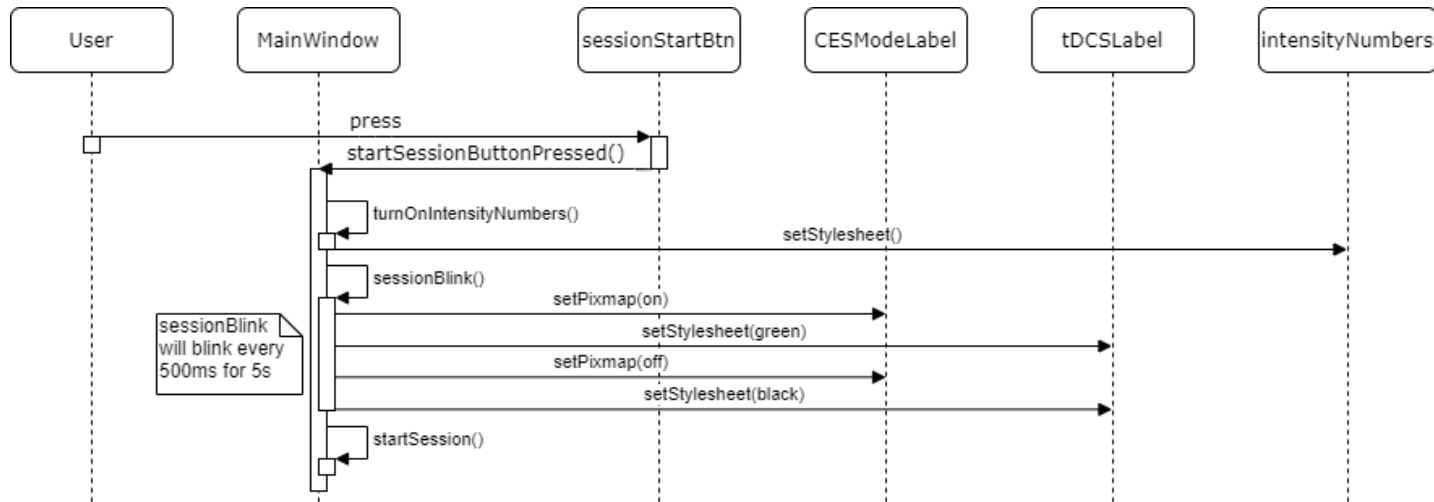


Adjusting Intensity

To adjust intensity, the user presses the INT buttons while a session is running to either increase the intensity or decrease the intensity

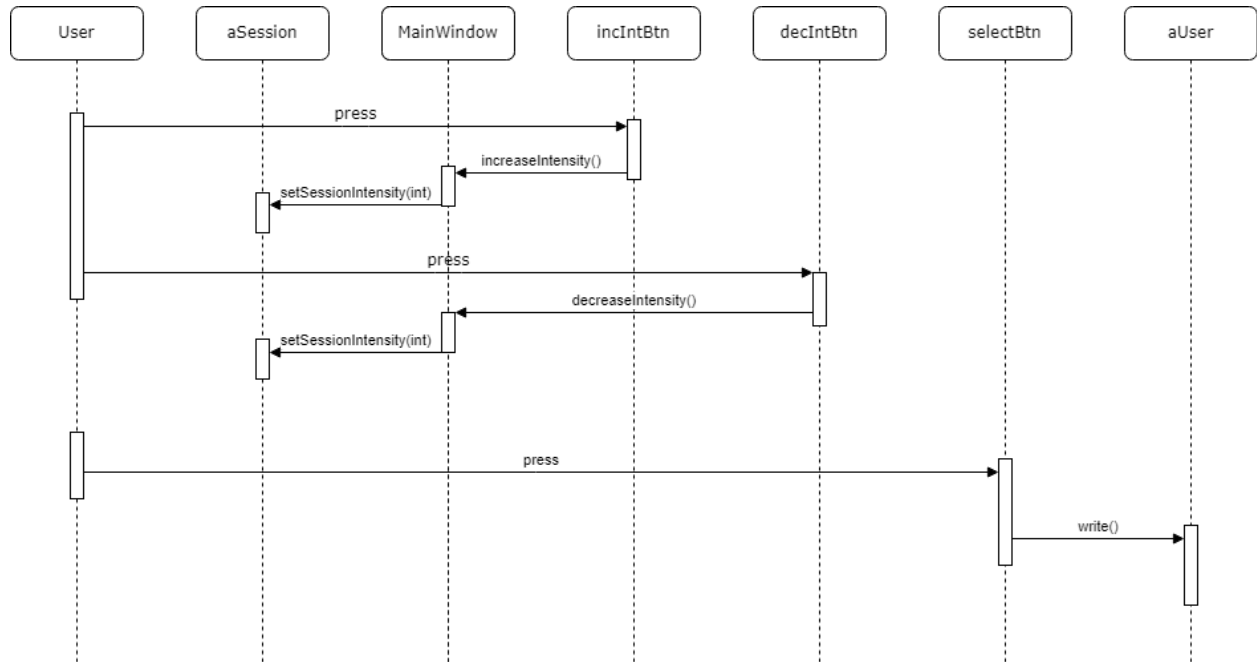

Switching users


To switch between the currently active user, press the switch user button located at the top right of the device, the currently selected user is highlighted beside the button.

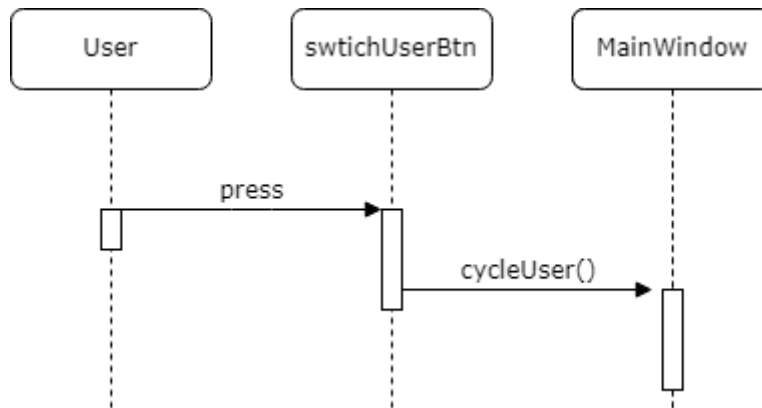

Recording a session

The user is able to record a session and save it for future use, to do so, the user first presses the record session button, which sets the machine into recording mode. Then the user can adjust the intensity, session type, and duration.

After the user has set their preferred settings, they can save it by pressing the save button.

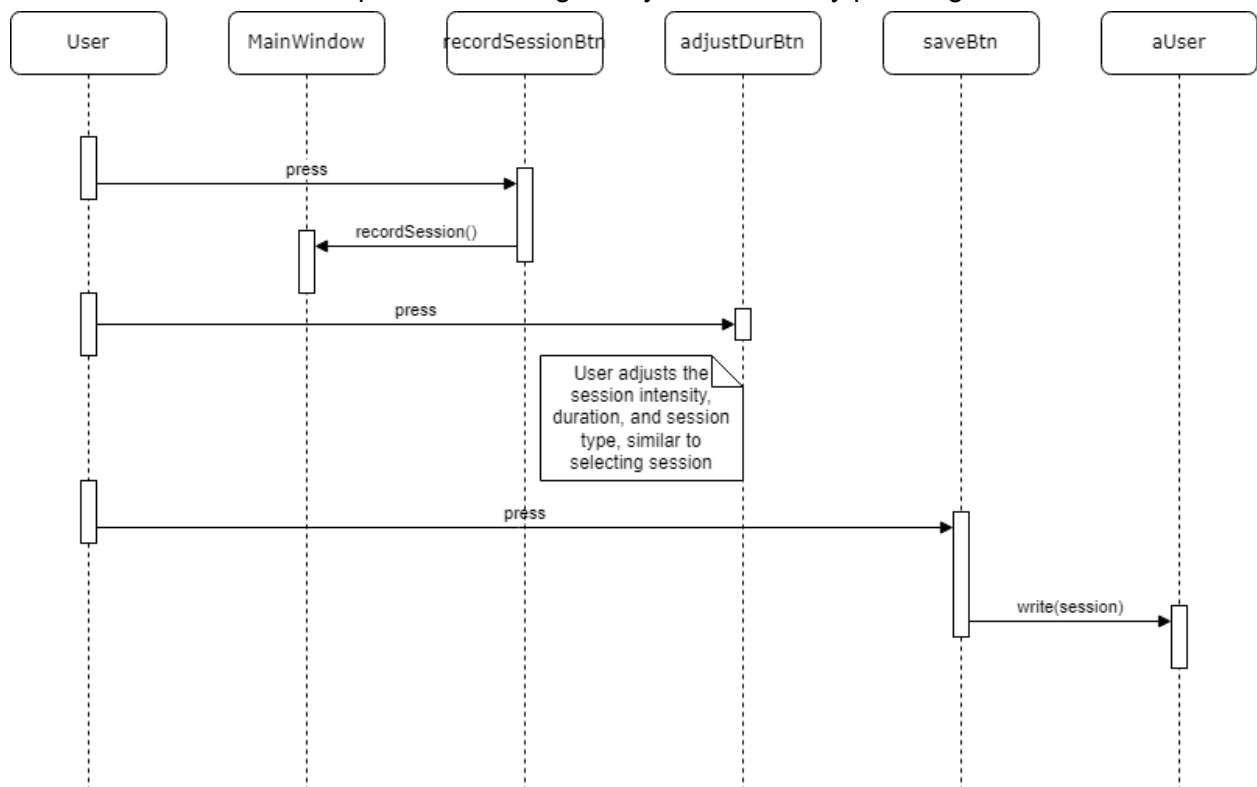

Replaying Session

To replay a session, the user selects a saved session from the session list which contains all saved sessions for the current user (limit of 4 saved sessions per user). The MainWindow then loads the selected session into memory. The user can then press the start session button to start their saved session.

